



GlyCLICK®

Azide Activation 250 µg

STORE AT

+4-8°C



FOR RESEARCH USE ONLY

Instructions for Use

GlyCLICK® Azide Activation 250 µg (L1-AZ1-025)
Process 250 µg IgG



DOWNLOAD INSTRUCTIONS FOR USE
www.genovis.com/ifu-L1-AZ1-025

Site-specific Conjugation of IgG with Azide-alkyne Click Chemistry

GlyCLICK is a site-specific conjugation technology for IgG using Fc glycan remodeling and click chemistry. The technology generates stable and homogenous antibody conjugates from several species and subclasses. Fc N-glycan remodeling by deglycosylation of the antibody allows for site-specific conjugation using robust click chemistry, resulting in a degree of labeling (DOL) or drug-antibody ratio (DAR) of 2.

GlyCLICK Azide Activation is available for azide activation of 250 µg, 2 mg or 10 mg IgG, for site-specific custom conjugation using an alkyne-modified label¹ of choice. All steps are performed under physiological conditions, thus maintaining the quality of the antibody. The site-specific conjugation on the Fc domain preserves the affinity of the antigen-binding sites.

GlyCLICK Azide Activation 250 µg contains all reagents needed to azide-activate 250 µg IgG. The azide activation is performed in two steps:

1. **Deglycosylation:** GlycINATOR Immobilized hydrolyzes the N-glycans on the Fc-part of the IgG to the inner GlcNAc.
2. **Azide Activation:** Azide attachment on the GlcNAc using GalT (Y289L)* and UDP-GalNAz*.

YOU MIGHT ALSO BE INTERESTED IN

GlycINATOR™ Immobilized

Immobilized enzyme for deglycosylation of IgG in spin columns

GlyCLICK® Biotin

Site-specific conjugation of IgG with biotin

GlyCLICK® Fluorophore


Site-specific conjugation of IgG with Alexa Fluor® 488, 555 or 647

GlyCLICK® DFO

Site-specific conjugation of IgG with DFO


1 For copper-free click chemistry, the label must be functionalized by a cyclooctyne.

* GalT (Y289L) and UDP-GalNAz are components of SiteClick™ and are provided under an intellectual property license from Life Technologies Corporation. The trademark SiteClick™ is the property of Life Technologies Corporation.

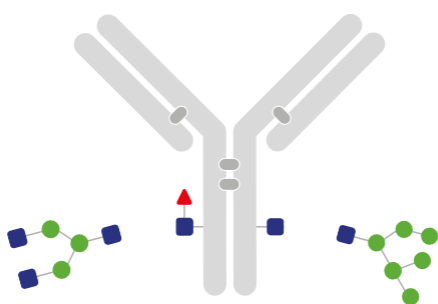

1. Deglycosylation

GlycINATOR Immobilized

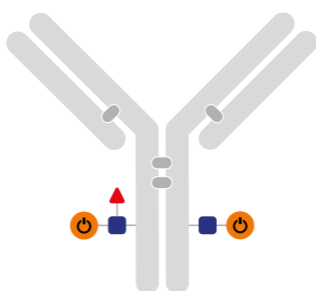

2. Azide Activation

GalT + UDP-GalNAz

Figure 1. Schematic overview of the GlyCLICK technology for azide activation.

Preparations

Important Information

Before you begin, briefly centrifuge tubes. Always wear suitable laboratory protective clothing and gloves when handling the reagents. **Keep in mind:** Sodium azide must be avoided throughout the protocol.

- Let the GlycINATOR Immobilized column equilibrate to room temperature before use.
- Use lids and bottom caps during the incubation.
- Before centrifugation, remove the bottom cap and loosen the lid (do *not* remove the lid).
- If a chelating agent will be used as a label, it is important to use metal-free water (trace analysis grade) throughout the protocol. The antibody must not be in contact with glass or metal.

Additional Materials Required

- IgG in 1×TBS, pH 7.4, free of carrier proteins and/or azide. 250 µg IgG in a maximum volume of 200 µl. To adjust the antibody solution, please follow “Guidance for Concentration and Buffer Exchange”. 20×TBS, a desalting spin column (0.5 ml, 40K) for buffer exchange and a concentrator (0.5 ml, 50K) are provided for convenience.
- 1×TBS: 10 ml 1×TBS is prepared by adding 0.5 ml of 20×TBS to 9.5 ml of ddH₂O. Vortex briefly to mix.
- Centrifuge tubes: 1.5-2 ml and 15 ml.
- ddH₂O.

Guidance for Concentration and Buffer Exchange

It is advisable to start with more IgG than 250 µg if concentration and/or buffer exchange of the sample is needed prior to “1. Deglycosylation: Modification of the N-glycan on the Antibody Fc Domain”.

A. Concentration Step

This step **is required if:**

- The volume of the IgG is more than 200 µl.

If the sample volume is 150-200 µl but needs a buffer exchange (if it contains phosphate or azide), concentrate the sample to < 150 µl and then follow the steps in section “B. Buffer Exchange”.

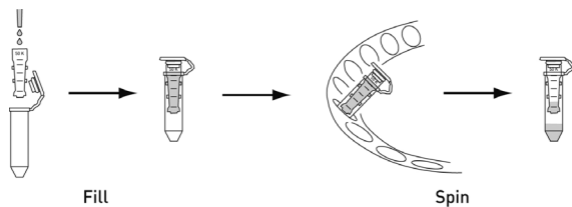


Figure 2. Antibody concentration step.

- A.1 Add 500 µl of ddH₂O to the antibody concentrator (0.5 ml, 50K) and cap the device as shown in Figure 2.
- A.2 Centrifuge at 5000 × g for 6 min. **Make sure that the cap strap and one membrane panel of the concentrator face the center of the rotor** (Fig. 2).
- A.3 Discard the flow-through.
- A.4 Add the IgG solution to the antibody concentrator.
- A.5 Centrifuge at 5000 × g for 2-6 min. **Make sure that the cap strap and one membrane panel of the concentrator face the center of the rotor** (Fig. 2).

Note: *If the IgG volume in the concentrator is more than 150 µl and the sample needs a buffer exchange, centrifuge for an additional 2 min at 5000 × g, or until the appropriate volume is achieved.*


- A.6 Invert the antibody concentrator into a new collection tube as shown in Figure 2.
- A.7 Centrifuge at 1000 × g for 3 min to collect the concentrated IgG. After collection, the volume of concentrated IgG should be approximately 100-150 µl in the collection tube.

B. Buffer Exchange

This step **is required if:**

- The IgG is in another buffer than 1x TBS, and/or
- The IgG is in a buffer containing azide

- B.1 Break off the bottom closure of the Desalting Spin column (0.5 ml, 40K). Loosen the lid (do not remove the lid).
- B.2 Place the column in a microcentrifuge tube and centrifuge at 700 × g for 1 min to remove the storage solution.
- B.3 Discard the flow-through and place the column in the microcentrifuge tube.
- B.4 Add 300 µl of 1× TBS on top of the resin. Centrifuge at 700 × g for 1 min and discard the flow-through.
- B.5 Perform step B.4 two additional times. The last time, centrifuge for 2 min.
- B.6 Blot the bottom of the column to remove excess liquid. Place the column in a new microcentrifuge tube.
- B.7 Apply the IgG solution (100-150 µl) on top of the resin.
- B.8 Centrifuge at 700 × g for 2 min and collect the flow-through containing the antibody in 1× TBS.

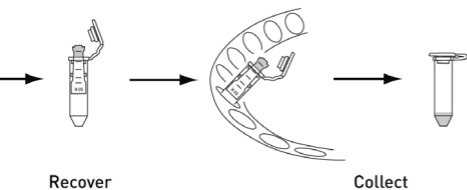

Recover

Collect

Site-specific Azide Activation of IgG

1. Deglycosylation: Modification of the N-glycan on the Antibody Fc Domain

Make sure that the antibody solution is 250 µg IgG in 100-200 µl of 1×TBS (pH 7.4, without azide).

Time required: 15 min hands-on, 30 min hands-off.

Materials from kit:

- 1× TBS (prepared from 20×TBS)
- Spin column with GlycINATOR Immobilized

- 1.1 Let the GlycINATOR Immobilized column equilibrate to room temperature before use. Break off the bottom cap of the GlycINATOR Immobilized column (save the cap) and place the column in a microcentrifuge tube. Loosen the lid.
- 1.2 Centrifuge at 200×g for 1 min to remove the storage solution. Discard the flow-through.
- 1.3 Place the column in the microcentrifuge tube.
- 1.4 Add 300 µl 1× TBS on top of the resin. Centrifuge at 200×g for 1 min and discard the flow-through.
- 1.5 Perform step 1.4 two additional times.
- 1.6 Insert the bottom cap.
- 1.7 Immediately add the antibody solution (100-200 µl) to the column.
- 1.8 Seal the column with the lid.
- 1.9 Fully suspend the media, mix by inversion and make sure there is a flow in the column.
- 1.10 Incubate the column with end-over-end mixing at room temperature for 30 min.
- 1.11 Remove the bottom cap and place the column in a new microcentrifuge tube. Loosen the lid.
- 1.12 Centrifuge at 1000×g for 1 min to collect the processed material.
- 1.13 Insert the bottom cap. Add 50 µl 1×TBS and seal the column with the lid.
- 1.14 Invert the column a couple of times.
- 1.15 Remove the bottom cap and place the column in a new microcentrifuge tube. Loosen the lid.
- 1.16 Centrifuge at 1000×g for 1 min to collect the processed material.
- 1.17 Repeat steps 1.13 to 1.16.
- 1.18 Pool the collected fractions, including the sample from step 1.12.

2. Azide Activation

Time required: 5 min hands-on, followed by overnight incubation.

Materials from kit:

- 1× TBS (prepared from 20× TBS)
- UDP-GalNAz
- GalT enzyme
- Buffer additive

- 2.1 Add 30 µl Buffer additive to the tube containing UDP-GalNAz.
- 2.2 Add the pooled fractions from step 1.18, and add 1× TBS to achieve a total volume of 375 µl.
- 2.3 Mix the solution by carefully pipetting up and down.
- 2.4 Add 25 µl of the GalT enzyme. The final reaction volume should be 400 µl. Mix the solution by carefully pipetting up and down. Seal the vial with the lid and wrap the lid with Parafilm® or similar.
- 2.5 Incubate overnight protected from light at 30°C.

3. Removal of Excess UDP-GalNAz

Time required: 1 h

Materials from kit:

- 1× TBS (prepared from 20× TBS)
- Antibody concentrator (2 ml, 50K)

- 3.1 Remove the conical collection tube from the antibody concentrator as shown in Figure 3.
- 3.2 Add 2 ml of 1× TBS to the antibody concentrator and centrifuge at 1200 × g for 10 min. **Make sure that one membrane panel of the concentrator faces the center of the rotor.**
- 3.3 Discard the flow-through.
- 3.4 Add 1.6 ml of 1× TBS and 400 µl of the azide-activated antibody from step 2.5 to the antibody concentrator (*Fig. 3*).
- 3.5 Centrifuge at 1200 × g for 6 min. **Make sure that one membrane panel of the concentrator faces the center of the rotor.**
- 3.6 Discard the flow-through.
- 3.7 Add 1× TBS to a total of 2 ml to the antibody concentrator.

Instructions for Use

3.8 Centrifuge at $1200\times g$ for 10 min. **Make sure that one membrane panel of the concentrator faces the center of the rotor.**

3.9 Discard the flow-through.

3.10 Perform steps 3.7-3.9 two additional times.

Note: *If the antibody volume in the concentrator is more than $\sim 200\mu\text{l}$, the volume in the concentrator can be reduced by additional centrifugation e.g. for an additional 5 min at $1200\times g$ or until the appropriate volume is achieved.*

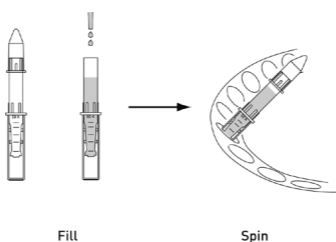


Figure 3. Use of the Antibody concentrator (2 ml, 50K).

3.11 Invert the antibody concentrator into the conical collection tube as shown in Figure 4.

3.12 Centrifuge at $1000\times g$ for 3 min to collect the concentrated azide-activated antibody.

3.13 Transfer the azide-activated antibody from the conical collection tube to a microcentrifuge tube.

3.14 At this stage, the azide-activated antibody can be stored at $2-8^{\circ}\text{C}$ protected from light for conjugation of a label at a later time.

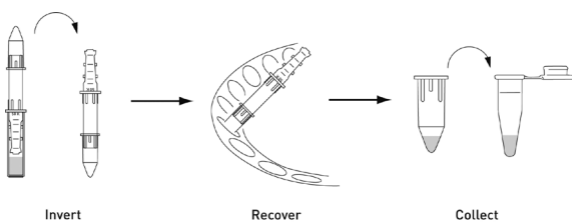


Figure 4. Collection of purified and concentrated azide-activated antibody.

CONTENT AND STORAGE

GlyCLICK Azide Activation 250 µg contains several components. The product box is shipped cold, and the components should be stored at +4-8°C upon arrival.

Do not freeze any of the components.

GlyCLICK Azide Activation is for R&D use only.

Table 1. Content and Storage Temperatures of GlyCLICK Components

Name	Amount	Store at
Antibody concentrator, 0.5ml, 50K	1 piece	4-25°C
Collection tube for antibody concentrator, 0.5ml, 50K	2 pieces	4-25°C
Desalting Spin column, 0.5ml, 40K	1 piece	4-8°C
GlycINATOR Immobilized, Microspin column	1 piece	4-8°C
UDP-GalNAz	1 vial solid	4-8°C Protect from light
20× TBS pH 7.4 (0.5M)	1 × 1.8 ml	4-8°C
Buffer additive	1 × 30 µl	4-8°C Protect from light
β-1,4-galactosyltransferase GalT (Y289L)	1 × 25 µl	4-8°C Protect from light
Antibody concentrator, 2ml, 50K	1 piece	4-25°C

USA & Canada

Genovis Inc.

10919 Technology Place Suite C, San Diego, CA 92127, USA

Phone: 1-855-782-0084 (toll free)

Fax: 1-858-524-3006

EMEA & Asia

Genovis AB

Box 4, SE-24421 Kävlinge, Sweden

Phone: +46 46 10 12 30

Fax: +46 46 12 80 20

support@genovis.com

www.genovis.com



All rights reserved. Genovis products may be covered by one or more patents, trademarks and copyrights owned or controlled by Genovis AB. For more information about commercial rights, please contact the Genovis team at licensing@genovis.com.

Genovis products are intended for research use only. They are not intended to be used for therapeutic or diagnostic purposes in humans or animals.

All goods and services are sold subject to Genovis' General Terms and Conditions of Sale.

SiteClick™ is provided under an intellectual property license from Life Technologies Corporation. The transfer of this product is conditioned on the buyer using the purchased product solely in research conducted by the buyer, excluding contract research or any fee for service research, and the buyer must not (1) use this product or its components for (a) diagnostic, therapeutic or prophylactic purposes; (b) testing, analysis or screening services, or information in return for compensation on a per-test basis; or (c) manufacturing or quality assurance or quality control, and/or (2) sell or transfer this product or its components for resale, whether or not resold for use in research. For information on purchasing a license to this product for purposes other than as described above, contact Life Technologies Corporation, 5791 Van Allen Way, Carlsbad, CA 92008 USA or outlicensing@lifetech.com.

The trademark SiteClick™ is the property of Life Technologies Corporation.

© Genovis AB